



Multichannel Sample Divider **Model 16:1**

Grain sample divider

Gehaka Multichannel Model 16:1



The ideal equipment for obtaining grain or seeds working samples without the manipulation of the operator.

The system works by gravity, without the need of any energy source.


It offers reliable and true results.

Reducing the working sample represents the initial batch, thus the Model 16:1 divider offers the certainty of a fair and impartial classification.

The channels have been engineered in such way that their dimensions are suitable to avoid obstructions, and the process of reducing the initial sample amount occurs in approximately 15 seconds.

The Gehaka Multichannel Sample Reducers can be used to receive grain, seeds and other granular products either from manually collected samples or from samples collected by pneumatic grain collectors systems.

Optional Accessories

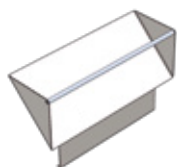


MULTIX Device - Allows making a second pass of the reduced product in the machine, which could increase or decrease the working sample.

(*)



Collector & Cone Dispenser - Enables the interconnection of the Gehaka Multichannel Homogenizer & Dividers with grain sampling's pneumatic collecting systems, allowing a continuous grain flow as a closed system.

(*)

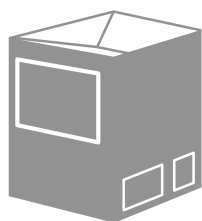

Available model: 16:1



Made in stainless steel,
aluminum and wood



Applications: grain and seed
receiving/shipping facilities,
laboratory rooms, universities,
grain processing facilities,
grain-warehousing facilities



Approx. hopper capacity:
Up to 15 kilograms of sample*.
*(when the optional Collector
& Cone Distributor is not used)
The final capacity will also depend
on the grain type and density

CONTACTS

gehaka.com.br/en

+55 (11) 2165-1100

vendas@gehaka.com.br

@gehakaoficial

/gehakaoficial

/gehakaoficial

/GehakaEquipamentos



Av. Duquesa de Goiás, 235 - Real Parque - CEP 05686-900 - São Paulo/SP - Brasil